



Yes, you can.[™]

PRESSURE CARE MATTRESSES



Welcome

Invacare is a global leader in care equipment, and we take pride in supplying quality products that meet the needs of both the user and carer.

Our Pressure Care Mattresses catalogue is designed to provide clear information to assist you in making a decision when choosing products that are essential in the provision of care in the home or within a facility.

It is not unusual for the needs of individuals to change and our comprehensive range of quality products are designed to provide solutions that meet many varying requirements of the individual.

Our global strength and our commitment to providing care, enables us to offer quality products and service that you can trust. Rely on us as a family member or as healthcare professionals providing support to you, and those you care for.

General Manager
Invacare Australia Pty Ltd



Visit our website

Log on to www.invacare.com.au regularly to keep updated on new product launches, and find out what events are on where you can view Invacare products.


Guide to Understanding Features

Invacare with its aim to making it simple for you to use our catalogue and products has provided a very simple key that gives you information at your fingertips.

Beds



Mattresses

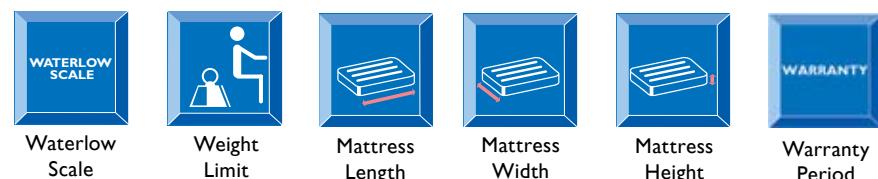


Symbol Guide

	Do not pierce or cut		Do not dry clean		Do not put near flame
	Wash at 80 °C		Do not bleach		User weight limit of 247 kg
	Cannot be tumble dried		Do not iron		Declaration of conformity

INVACARE®

Softform Premier Active**High/Very High Risk**

The Softform Premier Active retains properties as a high risk static mattress but should clinical judgement determine that the patient requires stepping up to a dynamic surface, then a discreet pump can be fitted to the air insert transforming the mattress into a dynamic surface and delivering additional levels of pressure relief.

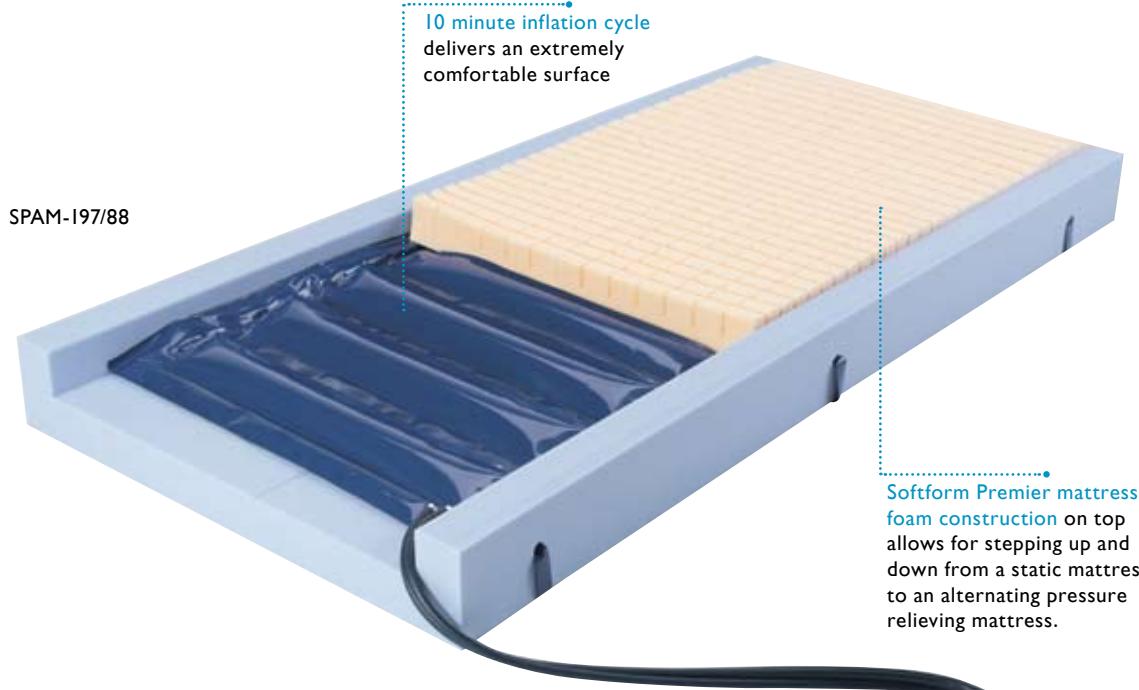



Features	Benefits
Air Insert contained within standard Softform Premier mattress foam construction	Allows stepping up and down between static and dynamic alternating systems with minimum patient disruption
Sensory software within the pump can gauge patient weight	Pressure can be maintained accordingly at an appropriate level. Fully alarmed in case of power failure/fault
Inflation of air cells alternates on a 10 minute cycle	Delivers an extremely comfortable mattress surface, while preventing pressure on any one area
Unique design of air channels	Will remain open in the event of constriction or rucking and maintain air pressure
Air system placed beneath foam insert	Greatly reduces risk of sharps puncture
Mattress can stay on the bed for a wide range of patients. Reduces requirement for introducing specialist dynamic systems	Reduces moving and handling issues; minimizes additional storage requirements
High risk static mattress, delivering exceptional pressure reduction and patient comfort	In the event of a power failure, the patient is still resting on a high risk pressure reducing surface, minimizing the potential of pressure damage occurring until power is restored
Mattress does not require turning	Reduces manual handling issues
Two-way stretch PU waterproof cover	Provides excellent recovery and durability, and helps reduce shear and friction forces
High frequency welded seams	Reduces the risk of fluid ingress
Full length flap conceals the zip	Reduces the possible fluid ingress through the zip
Toughened PU coated base	Prolongs the longevity of the mattress
Double coated PU sides	Reduces the risk of strikethrough, hydrolysis or delamination
Single piece of foam "U" core	Facilitates patient transfer by providing side wall support
Zipper for ease of removal	Unique four sided zip made from a moulded chain link plastic polymer that does not rust, whilst being extremely strong and durable
Inner surface light in colour	The white internal colour facilitates audits

Technical Specifications

	25		247.6kg		197cm		88cm		15.2cm		4 yr cover 8 yr foam
--	----	--	---------	--	-------	--	------	--	--------	--	-------------------------

INVACARE®

Softform Premier Glide

Very High Risk



Designed to conform to electric profiling beds, this mattress not only provides exceptional pressure reduction and patient comfort but its patented design allows the top surface supporting the patient to remain in position when profiling, significantly reducing the potential for shear and friction against patient tissue.

SPGM-197/88





Features	Benefits
Patented Glide system allows mattress to fully conform to the bed frame	Optimizes the properties of a profiling bed
Top surface remains under patient throughout profiling of an electric bed	Significantly reduces shear and friction forces, most particularly on the heel areas which are especially vulnerable to pressure ulcer development
Cover fabric conforms to deep contoured foam castellations on upper surface of the insert pad when mattress is profiled	Minimizes rucking of cover fabric on top surface, reducing potential for tissue damage
Does not require turning or rotating	Significantly reduces manual handling issues and helps to facilitate staff management resources
Two-way stretch PU waterproof cover	Provides excellent recovery and durability
High frequency welded seams	Reduces risk of fluid ingress and positively contributes to infection control
Full length flap conceals the zip	Reduces possible fluid ingress through the zip
Toughened PU coated base	Prolongs longevity of the mattress and reduces the risk of strikethrough, hydrolysis and delamination
Deep contoured insert pad for enhanced pressure area care	Maximise body contact area and optimizes pressure reduction
Single piece foam "U" core	Facilitates patient transfer by providing effective side-wall support
Does not impede CPR procedures	Clinically assessed to ensure this does not occur
Double-coated PU sides	Reduces potential for physical damage

Technical Specifications

	25		247.6kg		197cm		88cm		15.2cm		4 yr cover 8 yr foam
--	----	--	---------	--	-------	--	------	--	--------	--	-------------------------

INVACARE®


Softform Superior


SSMI97/88

Very High Risk

Mattress does not require turning reducing manual handling issues

The Superior mattress has been designed with tapered top surfaces, to improve pressure reduction in the high risk heel zone, and facilitate the easy use of the bed frame backrest.

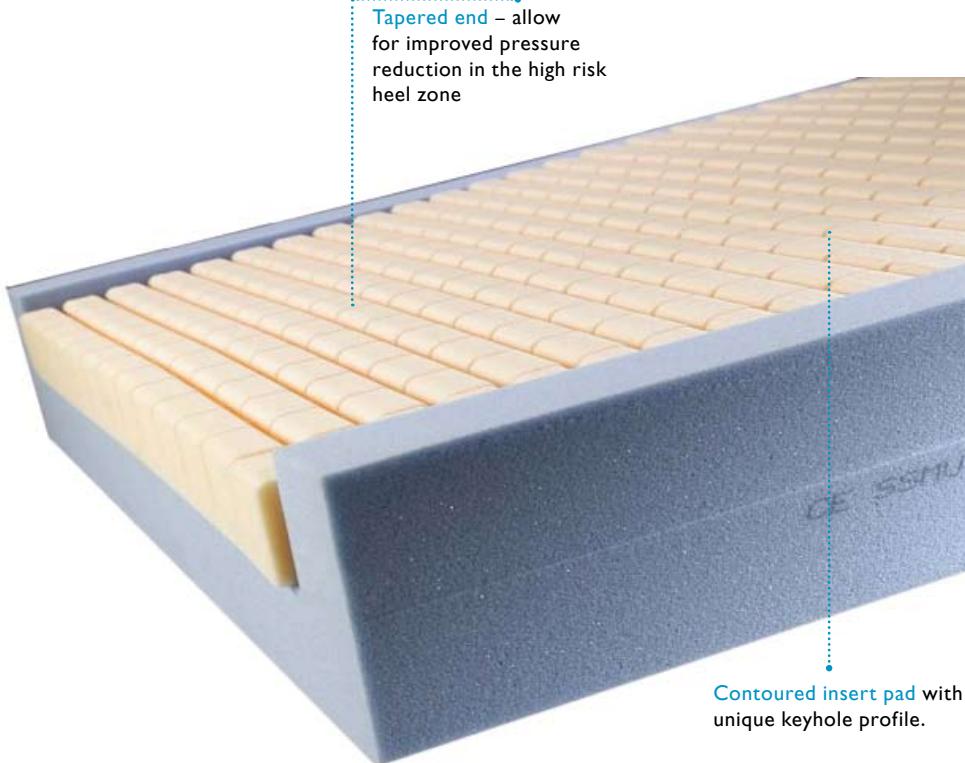

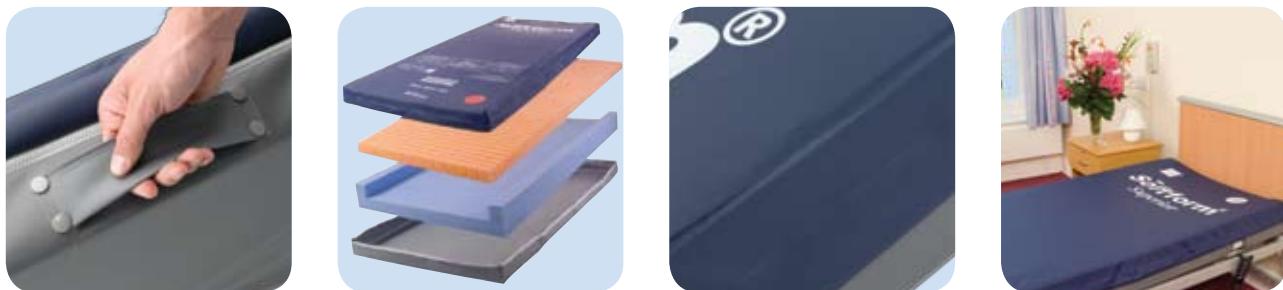

Features	Benefits
Tapering in heel zones	To improve pressure reduction in the high risk heel zone
Tapered top surfaces	Enable easy use of bed frame backrests
Cover fabric conforms to deep contoured foam castellations on upper surface of the insert pad when mattress is profiled	Minimizes rucking of cover fabric on top surface, reducing potential for tissue damage
Does not require turning or rotating	Significantly reduces manual handling issues and helps to facilitate staff management resources
Two-way stretch PU waterproof cover	Provides excellent recovery and durability
High frequency welded seams	Reduces possible fluid ingress through the zip
Toughened PU coated base	Prolongs longevity of the mattress and reduces the risk of strikethrough, hydrolysis and delamination
Deep contoured insert pad for enhanced pressure area care	Maximise body contact area and optimizes pressure reduction
Single piece foam "U" core	Facilitates patient transfer by providing effective side-wall support
Does not impede CPR procedures	Clinically assessed to ensure this does not occur
Double-coated PU sides	Reduces potential for physical damage
Top cover is fully detachable	Aids inspection and allows easy removal for washing. Also facilitates a cost benefit of allowing the replacement of individual parts rather than the whole mattress
Carrying handles	Improved transportation of the mattress (fitted to one side of the mattress only)

Technical Specifications

	25		247.6kg		197cm		88cm		15.2cm		4 yr cover 8 yr foam
--	----	--	---------	--	-------	--	------	--	--------	--	-------------------------

INVACARE®
Softform Premier


High Risk



Full length flap conceals
the zip reducing possible
fluid ingress

Equipped with a suite of features that ensures durability, sustained high degree of protection against tissue damage for the patient and an ingeniously designed insert pad, the Softform Premier mattress really does provide exceptional pressure reduction and comfort.

SPM-197/88

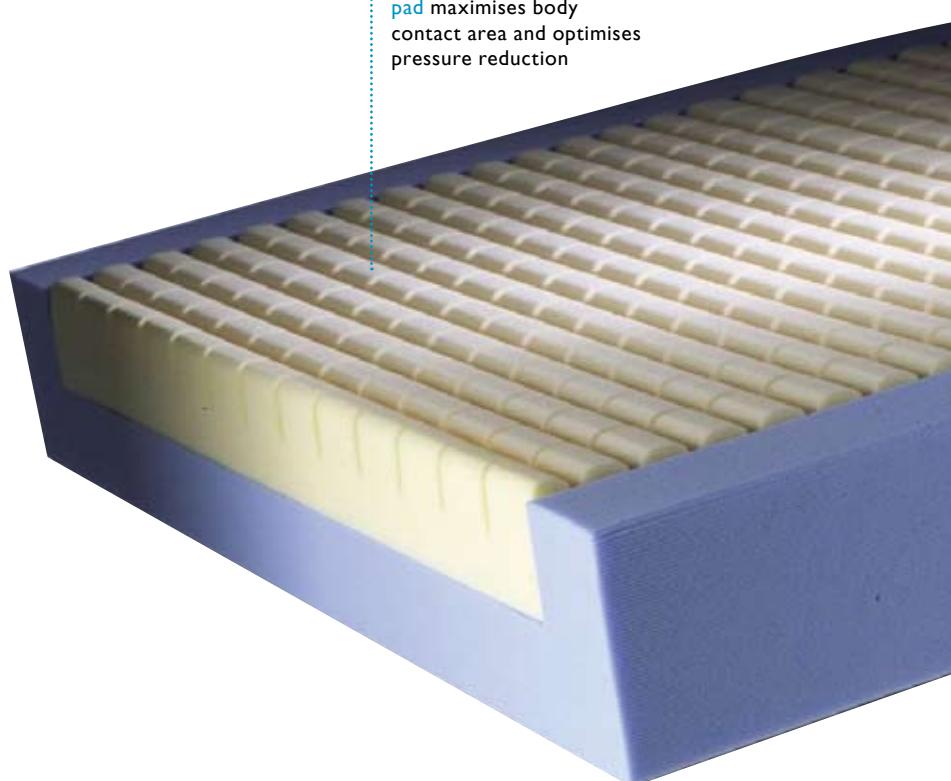

Deep contoured insert
pad maximises body
contact area and optimises
pressure reduction



Features	Benefits
Deep contoured insert pad for enhanced pressure care area	Maximises body contact area and optimises pressure reduction
Modular construction	Add an insert and pump at a later date to convert to a Softform Premier Active mattress
Cover fabric conforms to deep contoured foam castellations on the upper surface of the insert pad when mattress is profiled	Minimises rucking of cover fabric on the top surface, reducing potential for tissue damage
Suitable for High Risk patients / Very High Risk when used in conjunction with a profiling bed	Can significantly reduce pressure ulcer prevention
No need to turn over	Reduces manual handling issues
Single piece foam 'U' core	Facilitates patient transfer by providing effective side-wall support
Double-coated poly-urethane sides	Reduces the risk of strikethrough, hydrolysis or relamination
Two way stretch poly-urethane waterproof cover	Provides excellent recovery and durability and helps reduce shear and friction forces
High frequency welded seams and zip	Reduces risk of fluid ingress
Full length waterfall flap conceals the zip	Further risk reduction of fluid ingress (through zip)
Toughened PU coated base	Prolongs longevity of the mattress

Technical Specifications

	20		247.6kg		197cm		88cm		15.2cm		4 yr cover 8 yr foam
--	----	--	---------	--	-------	--	------	--	--------	--	-------------------------



HP100



HP50



OSW100

INVACARE®

Softform Accessories

INVACARE® Softform Heelpad

Designed to reduce pressure on the vulnerable heel area, the Softform Heel pad provides excellent immersion for high risk areas, the gel sacs allowing pressure displacement and even weight distribution.

The Softform Heel pads are available in double or single models. Both feature high frequency welded seams to prevent contamination, have a vapour permeable stretch cover with PU base and are comparatively lightweight (at 3.58kg and 1.5kg respectively) to enable safe and easy handling for carers.

INVACARE® Softform Odstock Wedge

Designed to reduce pressure under the heel, without increasing pressure on the thighs and calves the Odstock wedge helps support the patient from slipping down the bed, thus reducing shear forces in the sacral area.

By tilting the pelvis, it also reduces lumbar lordosis. This makes the Odstock wedge ideal after back surgery and improving the sleeping position of people with back problems.

Features	Benefits
Gel sac in Softform Heelpad	Providing excellent immersion for the high risk heel zone
High frequency welded seams	Protects products from contamination, ensuring product integrity
Vapour permeable stretch cover with PU base	High quality materials promote product longevity
Products designed to be used in conjunction to provide improved support to the lower limbs	Reduces the requirement for patient handling

INVACARE® Mattress Range


Invacares' range of mattresses are the intelligent choice when looking for a cost-effective mattress replacement option. Choose from a range that offers you safety and comfort for those at low to medium risk.



Comfort Economy

INVACARE® Comfort Economy Low Risk

The Comfort Economy is the economical option for those at Low Risk, a bi-layer convoluted foam mattress with an ultra thick base and low shear cover.

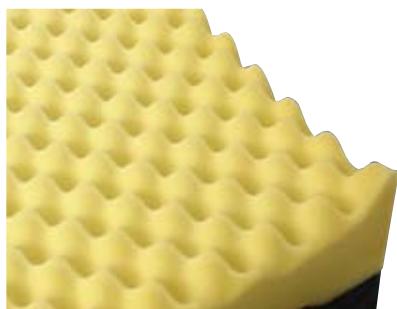

Comfort Plus

INVACARE® Comfort Plus Low Risk

The Comfort Plus is an economical and effective mattress featuring a convoluted top layer of foam on a bottom layer of solid foam and a low shear cover. It is ideal for low risk patients, reducing pressure and increasing comfort.

The base foam is an important part of the mattress structure. It prevents heavy patients from bottoming out, but provides a firm enough layer that cardio pulmonary resuscitation (CPR) can be carried out while the patient is on the mattress.



Foam Rest

INVACARE® Foam Rest **Low to Medium Risk**

The Foam Rest mattress has 2 layers of foam in a 2-way stretch cover. The top layer of foam features independent cubed cells which produce low shearing forces which allow for a more even distribution of pressure.



Low Resist

INVACARE® Low Resist **Low to Medium Risk**

The Low Resist pressure care mattress can be used for low to medium risk patients and is ideal for standard hospital or aged care beds. It provides excellent pressure care management and comfort to the user.



Safety

INVACARE® Safety **Low to Medium Risk**

The Safety range is an effective, multi-layered mattress available in low or low to medium risk models. Featuring solid foam "wings" on each side, a solid foam base and centre of either convoluted or cubed medium density foam in a low shear cover, the Safety has been designed to minimise the risk of falls by providing extra safety and support to the patient.*

* Affects Risk Level

INVACARE® Softform Mattress

Technical Specification

Product	Warranty	Fire Teasting	Grade Ref & Colour	Nominal Density Range (Kg/M)	Nominal Hardness Range (N)	Fatigue Class	Weight of Product (Kg)
Glide	4 Years Cover 8 Years Foam	BS 7177: Crib 5	Gr. RX 39/120 PINK Gr. RX 39/200 BLUE	38 - 40 38 - 40	105 - 135 180 - 200	Very Severe Very Severe	13
Superior	4 Years Cover 8 Years Foam	BS 7177: Crib 5	Gr. RX 39/120 PINK Gr. RX 39/200 BLUE	38 - 40 38 - 40	105 - 135 180 - 200	Very Severe Very Severe	13
Premier	4 Years Cover 8 Years Foam	BS 7177: Crib 5	Gr. RX 39/120 PINK Gr. RX 39/200 BLUE	38 - 40 38 - 40	105 - 135 180 - 200	Very Severe Very Severe	13
Premier Visco	4 Years*	BS 5852: Crib 5	Gr. RX 39/120 PINK Gr. RX 39/200 BLUE Gr. VC5580 WHITE	38 - 40 38 - 40 55	105 - 135 180 - 200 80	Very Severe N/A	13
Original	4 Years Cover 8 Years Foam	BS 7177: Crib 5	Gr. RX 39/120 PINK Gr. RX 39/200 BLUE	38 - 40 38 - 40	105 - 135 180 - 200	Very Severe Very Severe	12,5
Excel	2 Years	BS 7177: Crib 5	Gr. RX 39/120 PINK Gr. RX 39/200 BLUE	35 - 37 38 - 40	105 - 135 180 - 200	Very Severe Very Severe	10
A&E Trolley Mattress	1 Year	BS 5852: Crib 5	Gr. RX 36/125 BEIGE Gr. RX 39/200 BLUE	35 - 37 38 - 40	110 - 140 180 - 200	Very Severe Very Severe	5,5
Birthing Mattress	1 Year Cover 4 Years Foam	BS 7177: Crib 5	Gr. RX 39/200 BLUE	38 - 40	180 - 200	Very Severe	15
Newborn	1 Year	BS EN 597-1, 2	Gr. RN 26/100 WHITE	25 - 27	85 - 110	Average	0,5
Incubator Pad	1 Year	BS EN 597-1, 2	Gr. RN 26/100 WHITE	25 - 27	85 - 110	Average	0,5
Cot Mattress	1 Year	BS EN 597-1, 2	Gr. RN 26/100 WHITE	25 - 27	85 - 110	Average	2,7
Heel Pad	1 Year	BS 5852: Crib 5	Gr. RX 33/125 WHITE	32 - 34	110 - 140	Severe	3
Odstock Wedge	1 Year	BS 5852: Crib 5	Gr. RX 33/125 WHITE	32 - 34	110 - 140	Severe	2,5

* Standard Configuration

INVACARE® COMPANY PROFILE

INVACARE CORPORATION traces its roots back to 1885 when the Worthington Company began producing a line of vehicles designed for people with physical disabilities. Through the years, the complexion of the company has changed several times and the product range has expanded dramatically.

Today, Invacare is the global leader in the \$6 billion market for home medical products, manufacturing and distributing the world's broadest product offering to approximately 15,000 independent, home medical equipment providers.

Since establishment in Australia, Invacare has grown both organically and through the acquisition of several trusted and well known Australian brands including: Roller Chair; Pro-Med; Decpac; Australian Healthcare Equipment (Murphy Furniture and Hendicare Hospital Equipment); Morris Surgical; Home Health Equipment and Healthcare Equipment WA.

The acquisition of these Australian organisations has created a strong local business that has a broad product portfolio to support institutional and home healthcare.

Invacare Australia's current product range includes electric beds, furniture, static and active pressure care mattresses, patient lifters and slings, oxygen concentrators, manual and power wheelchairs, lightweight portable ramps, walking, toilet and bathroom aids.



Yes, you can.™

Invacare New Zealand
4 Westfield Place (PO Box 62-124) Mt Wellington Auckland, New Zealand
Freephone: 0508 468 222 Freefax: 0508 807 788
Email: sales@invacare.co.nz Web: www.invacare.co.nz

Invacare Australia Pty Ltd
1 Lenton Place (PO Box 5002) North Rocks NSW 2151, Australia
Freephone: 1800 460 460 Fax: 02 8839 5353
Email: sales@invacare.com.au Web: www.invacare.com.au

Invacare Asia Ltd
1 Lenton Place (PO Box 5002) North Rocks NSW 2151, Australia
Phone: +61 2 8839 5330 Fax: +61 2 8839 5331
Email: asiasales@invacare.com Web: www.invacare.com